


# InSite™

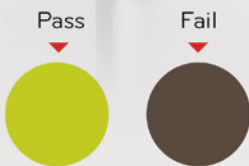


## Rapid Environmental Listeria Test

InSite is a quick and easy test for the detection of *Listeria* spp. in environments such as food processing facilities. The test contains a liquid media formulated with antibiotics, growth enhancers, and color changing compounds specific to the growth of *Listeria* species. A color change from yellow to light brown/black is considered presumptive positive. Applicability of InSite is limited for selected common *Listeria* spp (*Listeria monocytogenes*, *L. innocua*, *L. ivanovii*, and *L. welshmeri*) on selected common surface types.

**Key Features**

- Easy-to-use and no set up time
- Presumptive positive results are available in as little as 30 hours at levels as low as 1-10 CFU/mL
- Negative results after 48 hours
- Detects heat injured *L. monocytogenes* at low levels (10–50 cfu / ml of broth), even in the presence of high numbers of competing organisms
- Studies have shown equivalent or superior detection of injured *L. monocytogenes* when compared with standard culture media
- 12 month shelf life from date of manufacture
- Affordability - InSite cost 20 to 30% less than other rapid tests and reduces labor costs
- Only equipment required is a 37°C + 1 incubator

**RESULTS IN 30 HOURS**





**Americas**  
 941 Avenida Acaso  
 Camarillo, CA 93012  
 USA  
 Tel: 805.388.8007  
 Fax: 805.388.5531  
 Email: [enquiries@hygiena.net](mailto:enquiries@hygiena.net)  
 Web: [www.hygiena.net](http://www.hygiena.net)

**Europe**  
 Unit 25 Wenta Business Centre  
 Colne Way Watford  
 Hertfordshire WD2 4ND UK  
 Tel: +44 (0)1923 818821  
 Fax: +44 (0)1923 81825  
 Email: [enquiries@hygiena.net](mailto:enquiries@hygiena.net)  
 Web: [www.hygiena.net](http://www.hygiena.net)

**ORDER # IL50**

(50 tests per box)

**ORDER # IL100**

(100 tests per box)