


# QD-Loop™

## Labor Saving Dilution Devices

QD-Loop, a rapid-dilution device, is an all-in-one sterile system for accurate and convenient volume dilutions. An accurate calibrated loop comes attached to a reservoir where the patented Snap Valve™ technology device holds 1,0 ml of Butterfields or Maximum Recovery Diluent (MRD) buffer. Simply break the Snap Valve after a sample is collected and immediately a dilution is obtained that can be poured onto dry media film or streaked onto an agar plate. QD-Loop comes in 1/10 µl, 1/100 µl or 1/1000 µl dilution formats. This all-in-one system allows the making of dilutions quicker and more convenient.



Key Features	
✔	All-in-one sample collection system
✔	Multiple calibrated dilution loops
✔	Aseptically filled and gamma irradiated
✔	12-month shelf life
Benefits	
✔	No need to prepare dilution buffers
✔	No need to pipette
✔	Reduces labor and lab equipment expenses

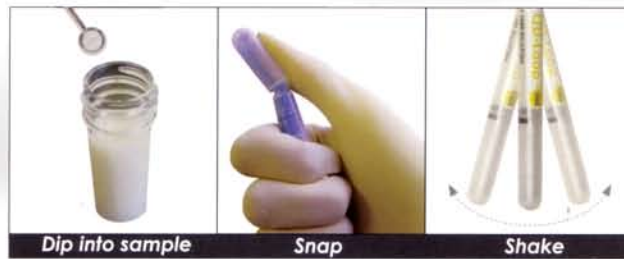

### Americas

941 Avenida Acaso  
 Camarillo, CA 93012  
 USA  
**Tel:** 805.388.8007  
**Fax:** 805.388.5531  
**Email:** enquiries@hygiena.net  
**Web:** www.hygiena.net

### Europe

Unit 11 Wenta Business Centre  
 Colne Way Watford  
 Hertfordshire WD2 4ND UK  
**Tel:** +44 (0)1923 818821  
**Fax:** +44 (0)1923 818825  
**Email:** enquiries@hygiena.net  
**Web:** www.hygiena.net

ORDER #	
<b>QD0010</b>	100µl sample pickup
<b>QD0100</b>	10µl sample pickup
<b>QD1000</b>	1µl sample pickup

(250 tests per box)