

Watershot™

Universal ATP Water Testing Device

Watershot is designed to help companies conducting water testing with a luminometer other than a Hygiena ATP system, reduce costs and achieve more accurate results. The 100µl dipper tip, combined with a unique liquid-stable reagent, provides more accurate sample collection and reproducible results than tests that use a swab bud. Watershot tests come in two different configurations:

Watershot Total - for detection of microbial and nonmicrobial ATP

Watershot Free - for detection of nonmicrobial ATP.



Microbial

Nonmicrobial

Key Features

- Simple, quick, and user-friendly
- Unique liquid-stable reagent provides unprecedented accuracy and reproducibility
- 12-month shelf life at refrigerated temperatures (2-8°C)
- 4-week shelf life at room temperatures (21°C)
- Tolerant to temperature abuse
- Collects 100µl of sample
- No need to change Pass/Fail setting - Watershot works exactly the same as the test it is replacing
- Costs up to 40% less than competitors ATP water tests

**COMPATIBLE WITH
THESE INSTRUMENTS**



ORDER #	Total	Free
3M / BioTrace	W-SPXL1333	WF-SPXL1333

(100 tests per box)

Related Products

Contamination Control Kit	-	CCK-3
Hold-Rite Cuvette Holder	-	HR001



Americas

941 Avenida Acaso
Camarillo, CA 93012
USA

Tel: 805.388.8007

Fax: 805.388.5531

Email: enquiries@hygiena.net

Web: www.hygiena.net

Europe

Unit 11 Wenta Business Centre
Colne Way Watford
Hertfordshire WD2 4ND UK

Tel: +44 (0)1923 818821

Fax: +44 (0)1923 818825

Email: enquiries@hygiena.net


Web: www.hygiena.net

Watershot Instructions

Watershot™ Total ATP Test with

Part No. W- SPXL1333

For Use With: 3M/Biotrace Unilite NG, Unilite XCEL, Unilite Luminometers

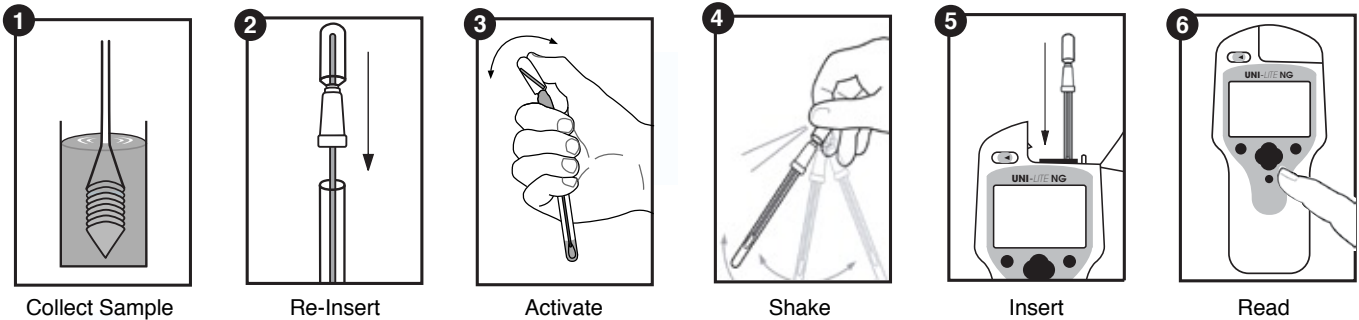

Watershot™ Free ATP Test with

Part No. WF- SPXL1333 | 100 Test

For Use With: 3M/Biotrace Unilite NG, Unilite XCEL, Unilite Luminometers

