

## QD-Loop™

### Labor-Saving Dilution Device


QD-Loop rapid dilution device is an all-in-one sterile system for accurate and convenient volumetric dilutions. An accurate calibrated loop comes attached to a reservoir where patented Snap Valve™ technology holds 1.0 mL of Butterfields or Maximum Recovery Diluent (MRD) buffer. Simply break the Snap Valve after a sample is collected and immediately a dilution is obtained that can be poured onto dry media film or streaked onto an agar plate. QD-Loop comes in 1/10, 1/100 or 1/1000 dilution formats. This all-in-one system allows the making of dilutions quicker and more convenient.

### Key Features

- All-in-one system
- Multiple calibrated dilution loops
- Aseptically filled and gamma irradiated
- 12-month shelf life

### Benefits

- No need to prepare dilution buffers
- No need to pipette
- Reduces labor and lab equipment expenses

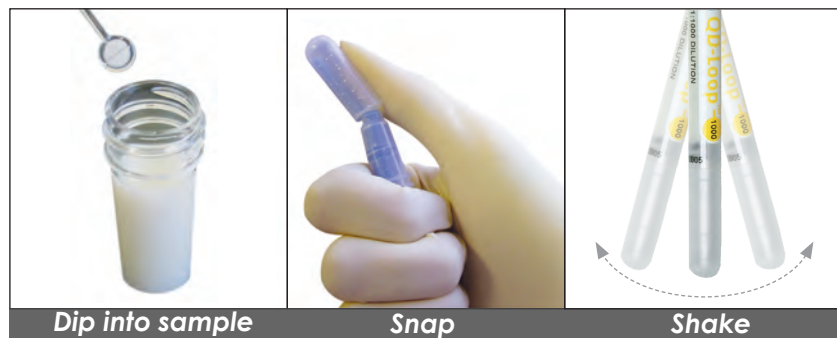

Order Numbers (250 per box):

QD0010 — 100µl Sample Pickup [1-10 dilution] \$176

QD0100 — 10µl Sample Pickup [1-100 dilution] \$176

QD1000 — 1µl Sample Pickup [1-1000 dilution] \$176



**hygiena**

Rapid Solutions for Food Safety

941 Avenida Acaso, Camarillo, CA 93012

Tel: 805.388.8007 ext. 300

email: [info@hygienaUSA.com](mailto:info@hygienaUSA.com)

web: [www.hygienaUSA.com](http://www.hygienaUSA.com)